

ZENITON 
Solution Kit


OSBC (OneSelf Bone Collector) Kit

Introduction & Feature

- Shaped to enable to collect the bone without any loss during drilling.
- Possible to do surgery with the quick & convenient drill classification, and selection due to color design per drill size.
- 4mm of stopper makes to drill safely.
- Possible length of drilling : 6mm, 8mm, 10mm
- Color division for identifying size
- Ideal design to remove bone



Direction for use

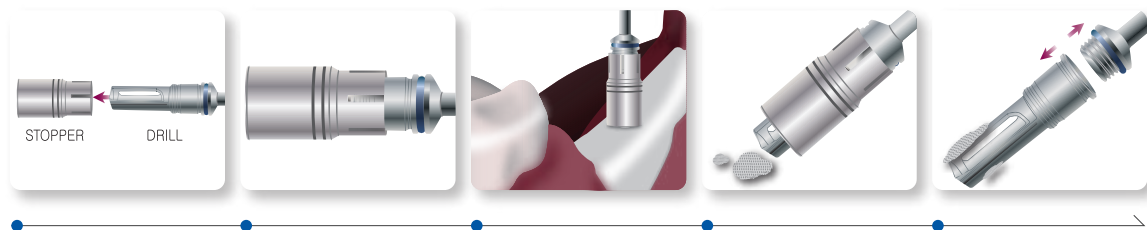

Clinical Case



OSBC (OneSelf Bone Collector)

Collector Drill



D	Ø3.5	Ø4.0	Ø4.5	Ø5.0
-	6501	6502	6503	6504

- Designed to collect the autogeneous bone during drilling, and the blade of drill can be separated from the shank.
 - Drill colors are different to distinguish the size of drill.
- [Recommended R.P.M. : 800 ~ 1,200 R.P.M.]

Stopper



L	4.0mm	6.0mm	8.0mm	10.0mm
-	6501S1	6501S2	6501S3	6501S4

- Up to 4mm, drilling can be done safely.
- Help with safe drill.

Bone Remover



-	BRM
---	------------

- Use for taking collected bone out after drilling.
- Use for removing drill from the shank.